



# Lympha-Flow

Sequential

Intermittent

Pneumatic

Compression



# Lympha-Flow (LF900)

Sequential Intermittent Pneumatic Compression

## Specifications

- + Power : 60VA
- + Time : 5~90min (unit : 5min)
- + Pressure : 10~180mmHg =-20mmHg (unit : 10mmHg)
- + Mode : Pre-therapy, Lymph, Distal Lymph, Soft lymph massage modes
- + Speed : 1~6 level
- + Dimension : 26(W)X18(D)X13(H)cm
- + Weight : 2kg (only main body)
- + Classification : Class II
- + Safety : IEC60601-1
- + EMC : IEC60601-1-12
- + Protection type : BF-type device

## Main functions

Operation

Continuous Use Setting

Manual Setting

Skip Setting

4 Protocol

Pre-Therapy Mode

Lymph Mode

Distal Lymph Mode

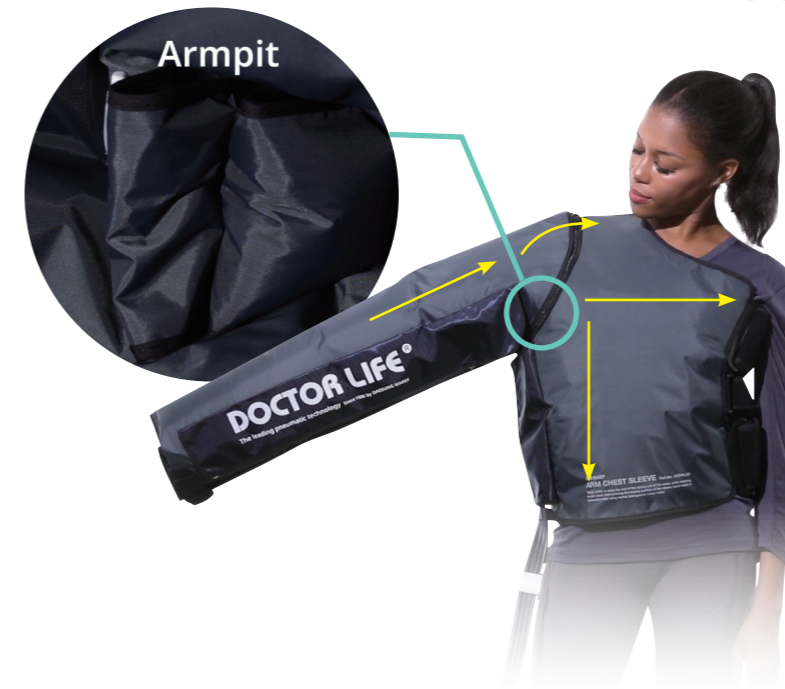
Soft Lymph Mode

Rental

Reset Function

Total Usage Time Record

Locking Program



## Upper Limb

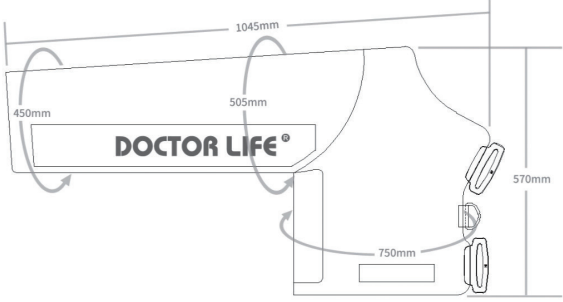
The ARM-CHEST garment has two specialized arm-pit chambers and nine overlapping chambers that help treat and move the stagnant lymph fluid in the correct physiological direct (into the upper body/ torso) after surgery of various orthopedic upper extremities such as breast cancer and breast plastic surgery.

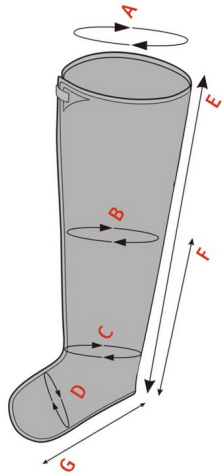
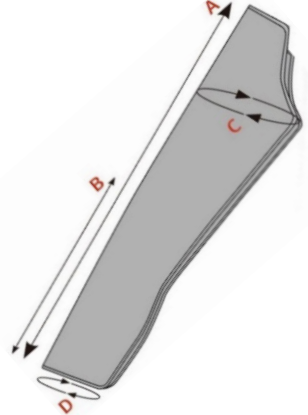
## Lower Limb

The lower limb sleeves such as leg sleeves and a center-body sleeve are designed to drain stagnant lymphatic fluid due to lymphatic damage after orthopedic surgery of various lower extremities. The center body sleeve can be more considered according to a doctor's prescription in case that the pelvic lymph node is damaged or abtained due to hysterectomy or ovarian cancer due to uterine leiomyoma or ovarian cancer, the stagnant lymphatic fluid occurred in the inguinal area.



# Overlapping Sleeves Dimensions

Ref, No	Descriptions	Dimension (mm)
4VDRL00	Arm Chest (Ac) Sleeve Medium	

Ref, No	Descriptions	A	B	C	D	E	F	G	Image		
S401	Full Leg Sleeve	S401 (Small)	760	640	555	435	600	200	285		
SW401		SW401 (Small Wide)	860	680	660	450	645	185	310		
SXW401		SXW401 (Small X-Wide)	990	890	730	460	600	210	310		
M401		M401 (Medium)	780	670	565	400	740	245	300		
MW401		MW401 (Medium Wide)	890	745	650	500	745	265	300		
MXW401		MXW401 (Medium X-Wide)	1020	880	750	545	700	240	320		
L401		L401 (Large)	790	680	520	440	860	290	300		
LW401		LW401 (Large Wide)	1000	720	560	4420	890	285	290		
M403		Half Leg Sleeve	M403 (Medium)	760	630	550	440	420	130		300
MW403			MW403 (Medium Wide)	810	675	660	580	490	165		300
L402	Arm Sleeve	M402 (Medium)	700	360	560	330	x	x	x		
L402		L402 (Large)	820	410	660	420	x	x	x		